

Breezeway® II

DELIVERY SHEATH WITH LOADER

Making Transeptal Crossing a Breeze-way!™

Breezeway® II


The Breezeway® II Delivery Sheath is intended for the introduction of diagnostic and therapeutic devices into the human vasculature, including intracardiac.



DELIVERY SHEATH TIP & DILATOR FEATURES


- Soft atraumatic tip
- Smooth tapered sheath tip dilator transition during peripheral and transseptal access
- Superior curve retention for reliable catheter placement
- Sideport holes on distal tip for aspirating and flushing the deliver sheath
- The dilator is constructed of a low friction and radiopaque material

LEFT ATRIUM ACCESS

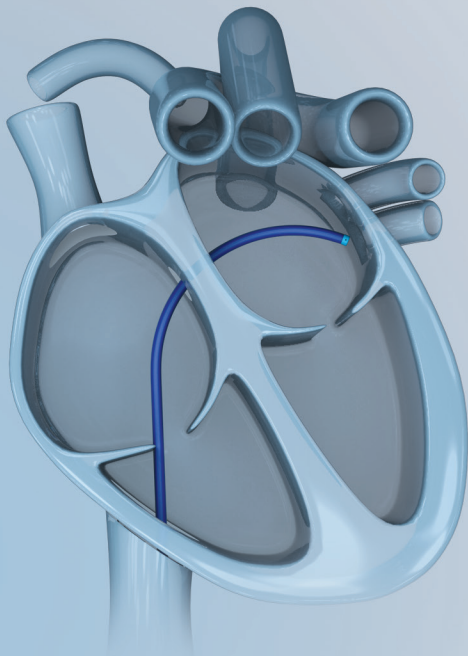

Dilator

Sheath

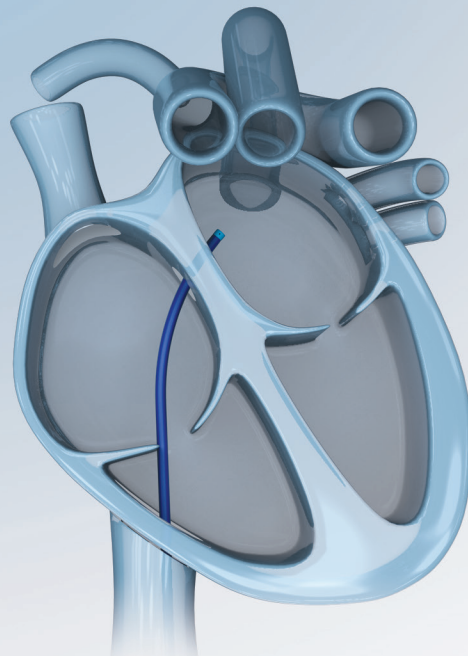


Dilator

Sheath

DELIVERY SHEATH SHAFT FEATURES

- Hydrophobic coated sheath
- Low friction nylon inner liner
- Clear view of the distal marker band under fluoroscopy for easy positioning

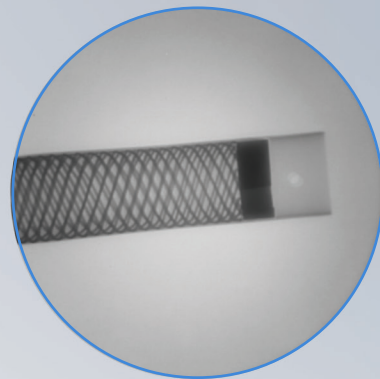

FLEXSTEER™ PROPRIETARY SHAFT DESIGN


- Reliable Flexsteer™ reinforced braid technology provides exceptional control, kink resistance, while maintaining optimal flexibility

Liner

Braid

Polymer

Coating

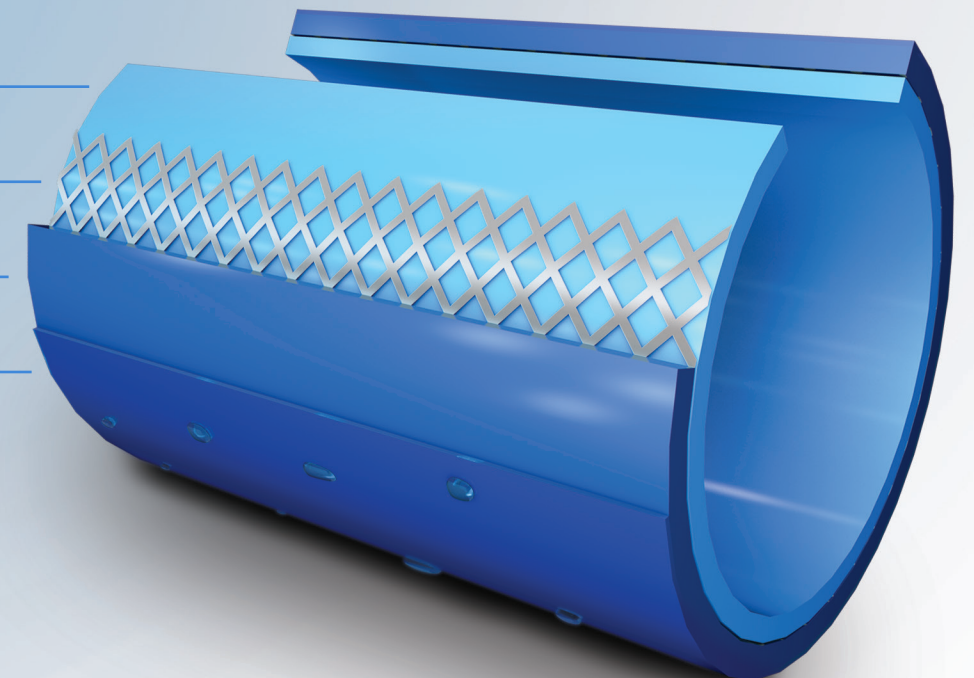

CONNECT & DELIVER



DELIVERY SHEATH HUB FEATURES


- Optimal SureSeal™ hemostatic valve technology to minimize air embolism and blood loss
- Available in two hub designs 6F - 14F

LOADERS




LOADER FEATURES

- Loader compatible hub connector, available for all French sizes
- Distal hub lock to attach the loader to the sheath
- Transparent material to allow visualization of air bubbles
- Low friction inner lumen
- Available in two hub designs 6F - 14F

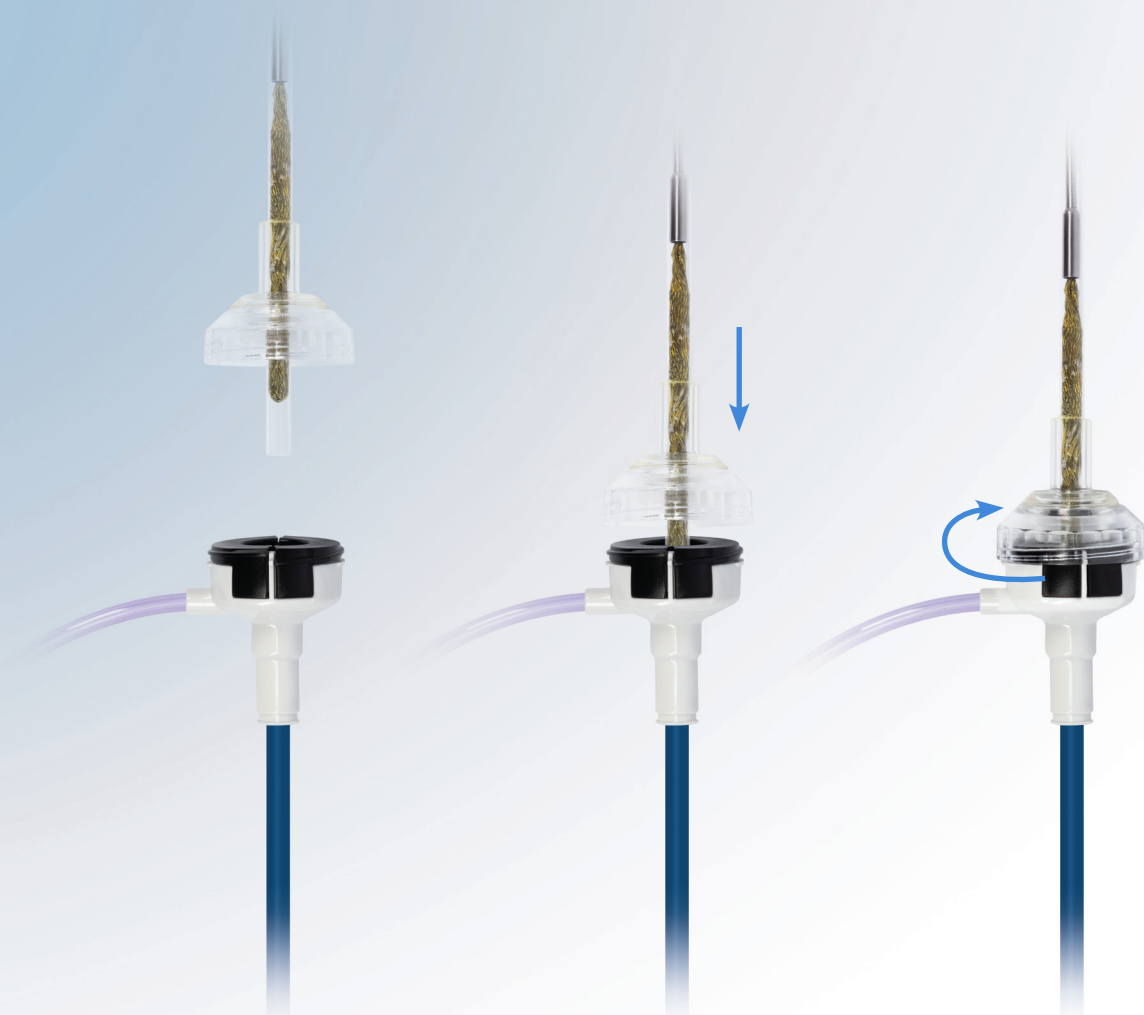


ORDER INFORMATION

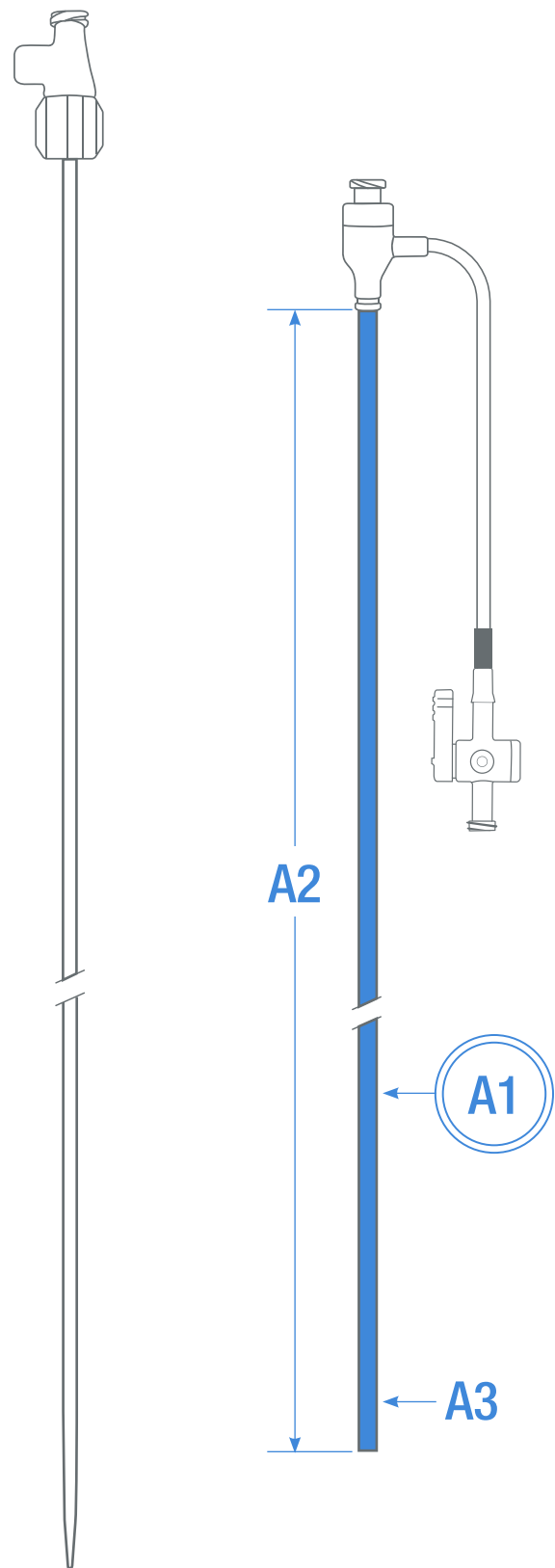


- Kit Includes:**
- (1) Breezeway® II Sheath
 - (1) Breezeway® II Dilator

- Plus Kit Includes:**
- (1) Breezeway® II Sheath
 - (1) Breezeway® II Dilator
 - (1) Loader

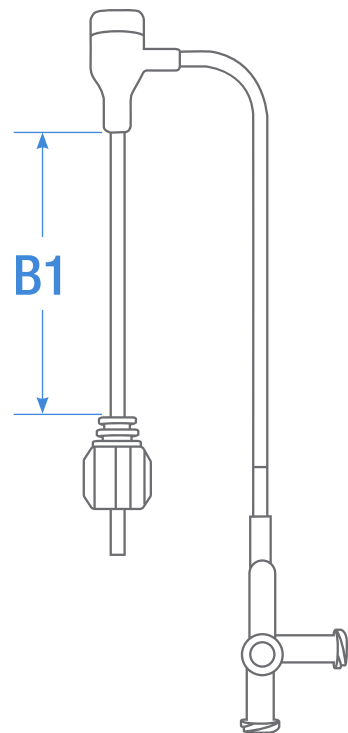

SHEATH DIAMETERS

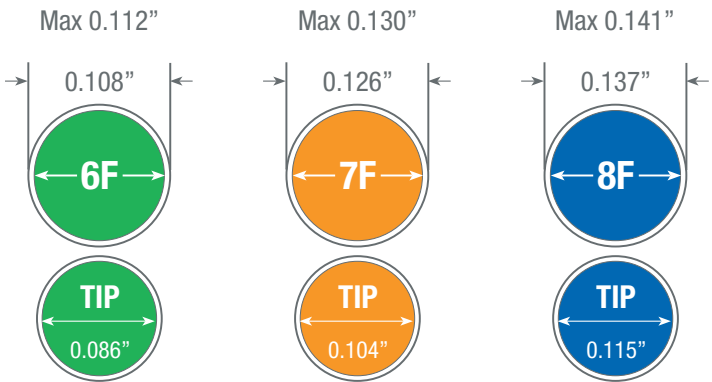

LOADER DIAMETERS




BREEZEWAY® II SPECIFICATIONS					
Model # Kit (One per box)	Model # Plus Kit (One per box)	A1 Sheath Inner Diameter	A2 Usable Sheath Length	A3 Sheath Tip Curve	B1 Usable Loader Length
BG060079045	BG060079045L	● 6 F	79 cm	45°	8.5 cm
BG060110045	BG060110045L	● 6 F	110 cm	45°	8.5 cm
BG070079045	BG070079045L	● 7 F	79 cm	45°	8.5 cm
BG070079180	BG070079180L	● 7 F	79 cm	180°	8.5 cm
BG070110045	BG070110045L	● 7 F	110 cm	45°	8.5 cm
BG080079045	BG080079045L	● 8 F	79 cm	45°	12.5 cm
BG080079180	BG080079180L	● 8 F	79 cm	180°	12.5 cm
BG080110045	BG080110045L	● 8 F	110 cm	45°	12.5 cm

*Delivery sheath supplied sterile. For single use only. Available bulk upon request.
Standard lengths are noted above. Other lengths and curves available upon request.*


Products not available in all regions.

ORDER INFORMATION

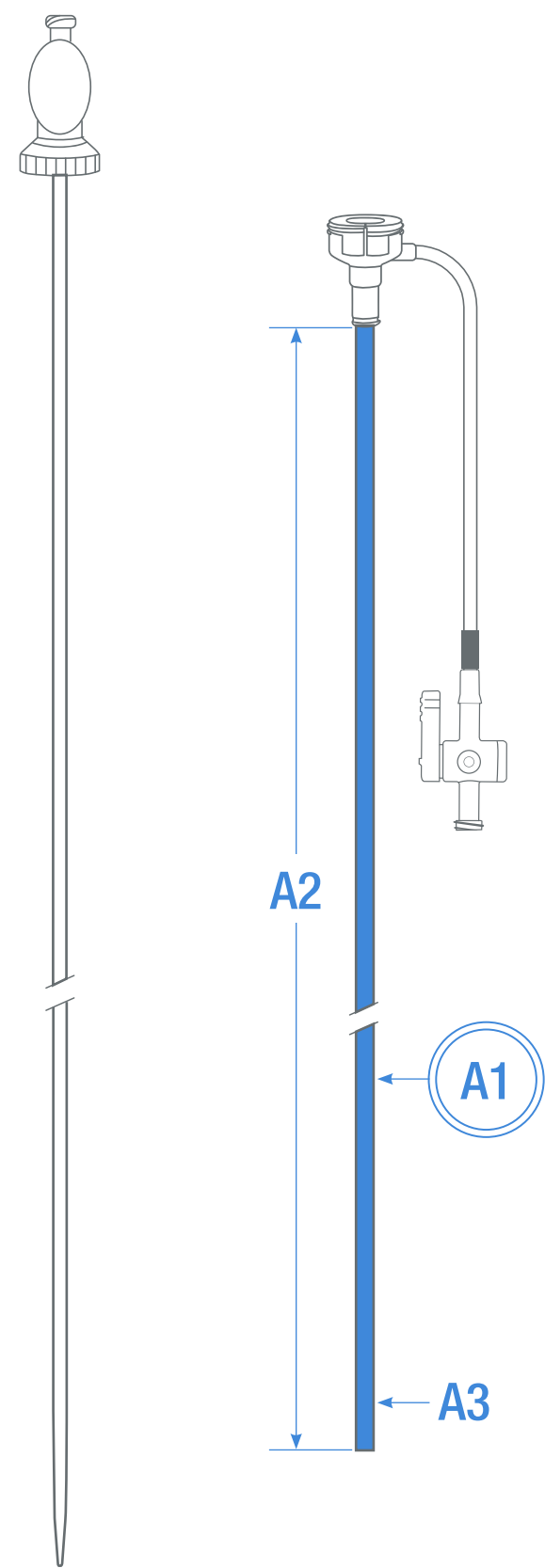


- Kit Includes:
- (1) Breezeway® II Sheath
 - (1) Breezeway® II Dilator

- Plus Kit Includes:
- (1) Breezeway® II Sheath
 - (1) Breezeway® II Dilator
 - (1) Loader

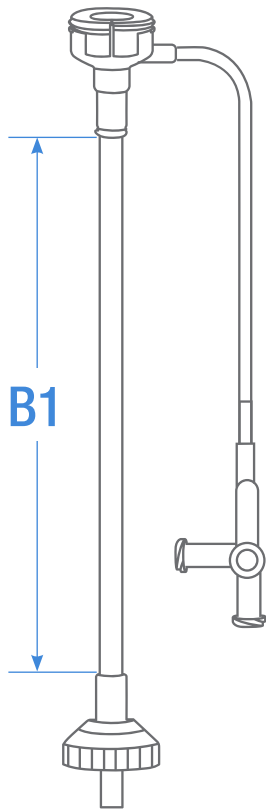

SHEATH DIAMETERS

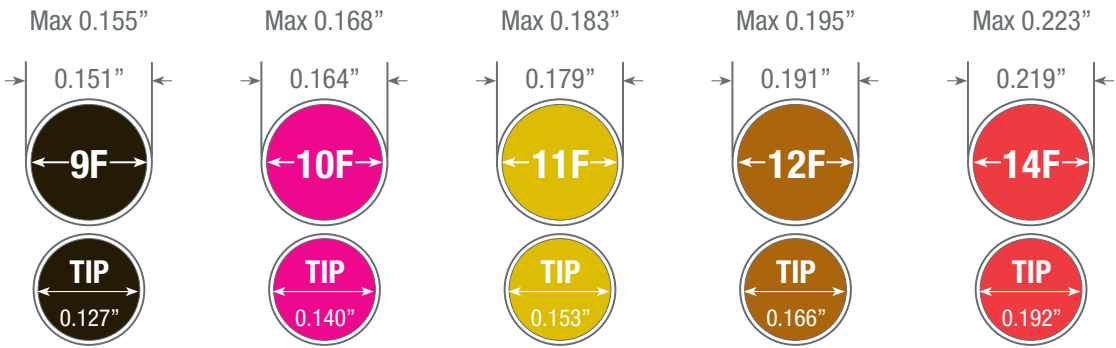

LOADER DIAMETERS



BREEZEWAY® II SPECIFICATIONS					
Model # Kit (One per box)	Model # Plus Kit (One per box)	A1 Sheath Inner Diameter	A2 Usable Sheath Length	A3 Sheath Tip Curve	B1 Usable Loader Length
BG090079045	BG090079045L	● 9 F	79 cm	45°	12.5 cm
BG090079180	BG090079180L	● 9 F	79 cm	180°	12.5 cm
BG090110045	BG090110045L	● 9 F	110 cm	45°	12.5 cm
BG100079045	BG100079045L	● 10 F	79 cm	45°	12.5 cm
BG100079180	BG100079180L	● 10 F	79 cm	180°	12.5 cm
BG100110045	BG100110045L	● 10 F	110 cm	45°	12.5 cm
BG110079045	BG110079045L	● 11 F	79 cm	45°	16.5 cm
BG110079180	BG110079180L	● 11 F	79 cm	180°	16.5 cm
BG120079045	BG120079045L	● 12 F	79 cm	45°	16.5 cm
BG140079045	BG140079045L	● 14 F	79 cm	45°	16.5 cm

Delivery sheath supplied sterile. For single use only. Available bulk upon request.
Standard lengths are noted above. Other lengths and curves available upon request.

Products not available in all regions.



OSCOR INC. | WORLD HEADQUARTERS

3816 DeSoto Blvd. Palm Harbor, FL 34683 – USA
Toll Free: 1-800-726-7267
P: +1 727-937-2511

OSCOR EUROPE GMBH

Telefon kostenlos: 0800 00-67267
P: +49 (0) 211 586 786-00

CONTACT US

www.oscor.com
sales@oscor.com

