

Baylis
MED TECH

TWIST

UNI-DIRECTIONAL
DEFLECTABLE
GUIDING SHEATH

DESTINO™


TWIST

A New Twist On Access Technology




Maximizing Sheath Control

- 180° deflectable distal sheath tip section allows smooth steering and navigation through the vascular system
- Clear view of the distal marker band under fluoroscopy for easy positioning

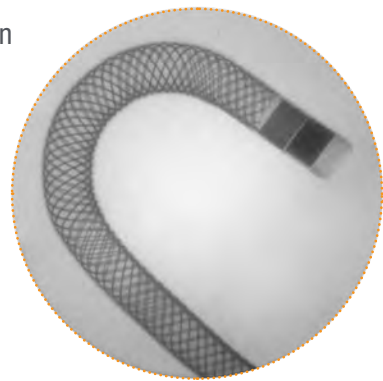

- PTFE Liner
- Flexsteer™ reinforced braid
- Outer wall extrusion
- Hydrophobic coated sheath

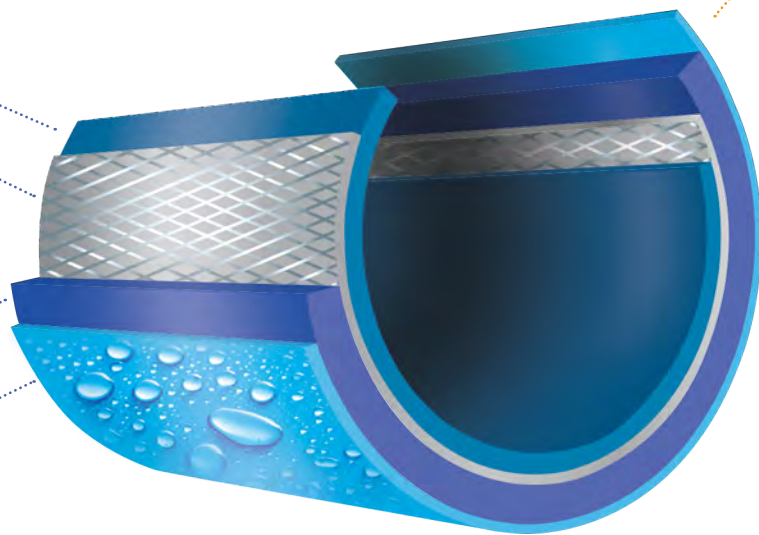

- Flexsteer™ reinforced braid technology provides improved torque control and kink resistance, while maintaining optimal flexibility even in the most tortuous anatomy


Minimizing Sheath Exchanges

YOUR ONE SHEATH SOLUTION

- The Destino Twist deflectable tip minimizes your need for multiple fixed curves and catheter exchanges




- SureSeal™ hemostatic valve technology minimizes air embolism and blood loss

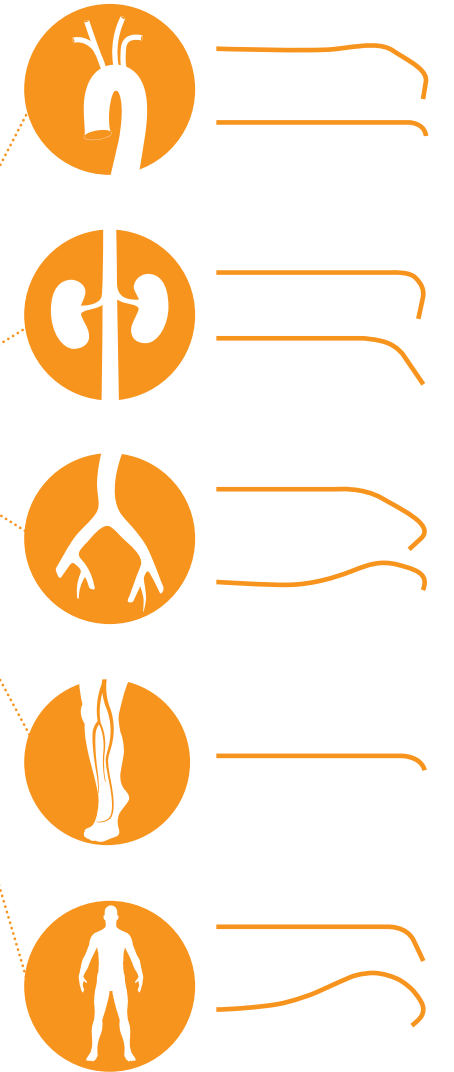

Maximizing Your Options

Diameter Options:

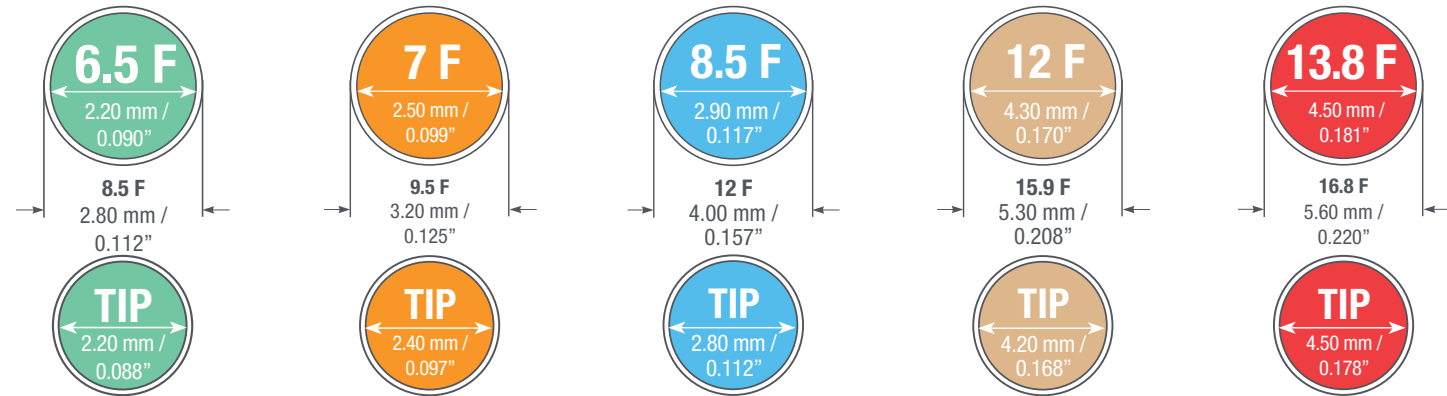

Curve Options:



- Clear view of the distal marker band under fluoroscopy
- Smooth tip to dilator transition
- Hydrophobic coated sheath
- Braided Flexsteer™ shaft technology
- Rotating collar
- Tip position indicator
- Reliable SureSeal™ hemostatic valve
- Dilator with French size and guidewire indicator

12 F / 13.8 F dilator cap Standard dilator cap


12 F / 13.8 F seal Standard seal



CAROTID ACCESS

A New Twist on Stability.

The Destino Twist simplifies access to challenging arch anatomy. The deflectable tip gives the user a stable platform and minimizes device exchanges.

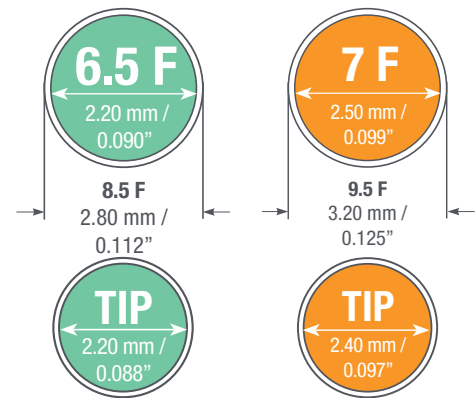


Model # Set <i>(One per box)</i>	Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend
STANDARD				
DST0659009	6.5 F	90 cm	108 cm	9 mm
DST0709009	7 F	90 cm	108 cm	9 mm

ILIAC/SFA ACCESS

A New Twist on Advanced Support.

The Destino Twist uni-directional sheath features a variety of curves, atraumatic tip, and hydrophobic coating which improves trackability. The Twist's ability to deflect, anchors and stabilizes the sheath across the bifurcation enhancing contralateral access.

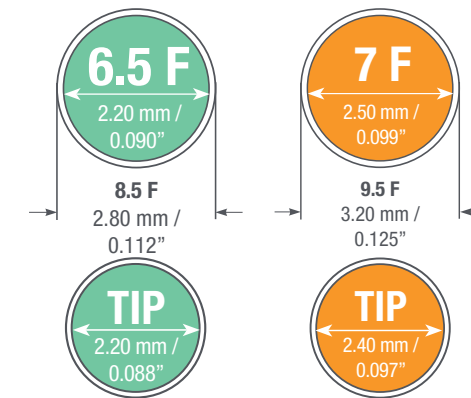

Model # Set <i>(One per box)</i>	Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend
STANDARD				
DST0654509	6.5 F	45 cm	63 cm	9 mm
DST0654517	6.5 F	45 cm	63 cm	17 mm
DST0655507D2*	6.5 F	55 cm	73 cm	7 mm
DST0704509	7 F	45 cm	63 cm	9 mm
DST0704517	7 F	45 cm	63 cm	17 mm


* Contains a second soft dilator with 80 cm usable length.





INFRAINGUINAL/ INFRAPOPLITEAL ACCESS

A New Twist on Deflection Control.

The Destino Twist offers stable platform for access below the knee. Used with a coaxial selective catheter, it simplifies access while minimizing device exchanges.

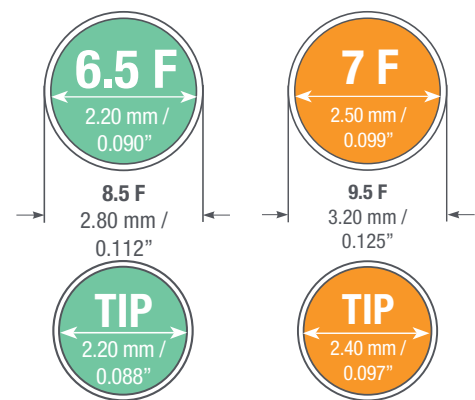

Model # Set <i>(One per box)</i>	Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend
STANDARD				
DST0654509	● 6.5 F	45 cm	63 cm	9 mm
DST0654517	● 6.5 F	45 cm	63 cm	17 mm
DST0704509	● 7 F	45 cm	63 cm	9 mm
DST0704517	● 7 F	45 cm	63 cm	17 mm

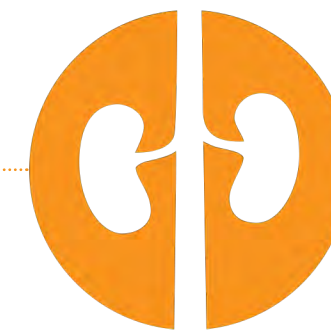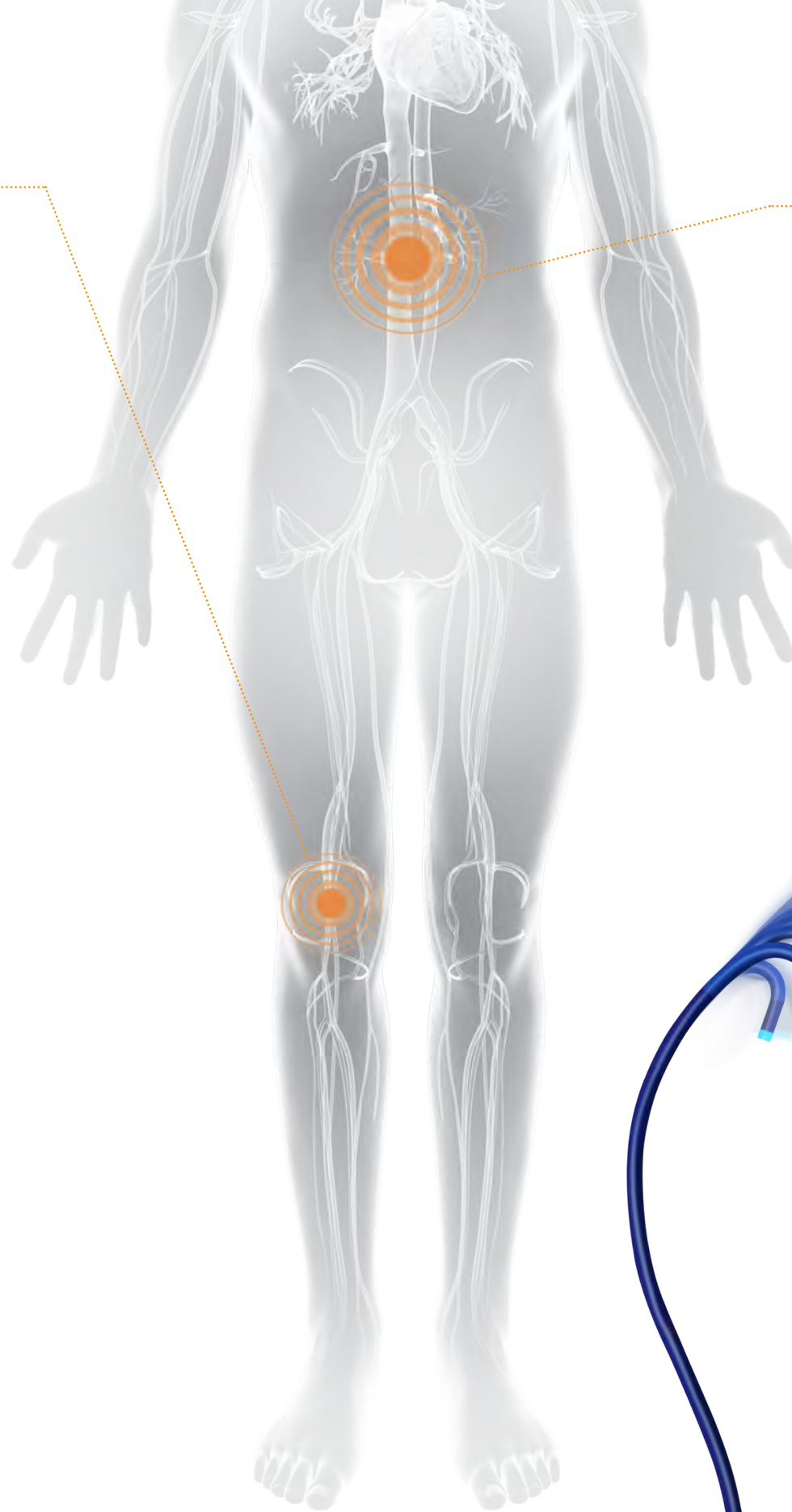


RENAL/CELIAC/SMA ACCESS

A New Twist on Accessing Challenging Anatomy.

The Destino Twist's deflectable tip simplifies access to difficult patient anatomies, which is especially useful when accessing arteries with inferior take offs. This enables the physician to access from the groin instead of using a brachial approach.

Model # Set <i>(One per box)</i>	Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend
STANDARD				
DST0654509	● 6.5 F	45 cm	63 cm	9 mm
DST0654517	● 6.5 F	45 cm	63 cm	17 mm
DST0655509	● 6.5 F	55 cm	73 cm	9 mm
DST0655517	● 6.5 F	55 cm	73 cm	17 mm
DST0704509	● 7 F	45 cm	63 cm	9 mm
DST0704517	● 7 F	45 cm	63 cm	17 mm
DST0705509	● 7 F	55 cm	73 cm	9 mm
DST0705517	● 7 F	55 cm	73 cm	17 mm
TD				
DST0655509TD	● 6.5 F	55 cm	73 cm	9 mm
DST0705509TD	● 7 F	55 cm	73 cm	9 mm


TD (*Tichà Damascelli*)
Curve for added
back support

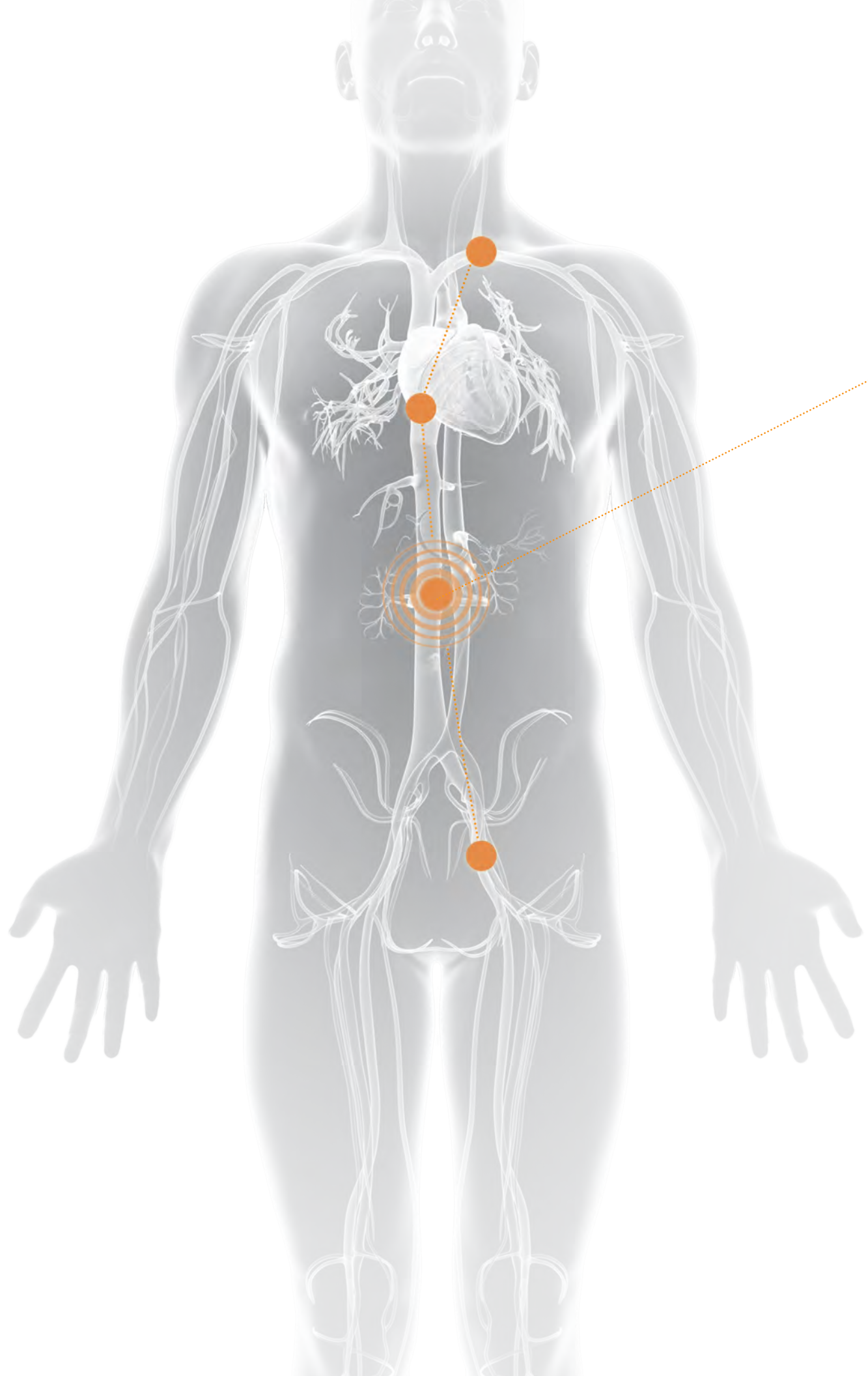

MULTIPURPOSE ACCESS

A New Twist on Adaptable Deflection.


The Destino Twist diverse product offering allows for optimal compatibility with patient anatomy adapting to all your needs.

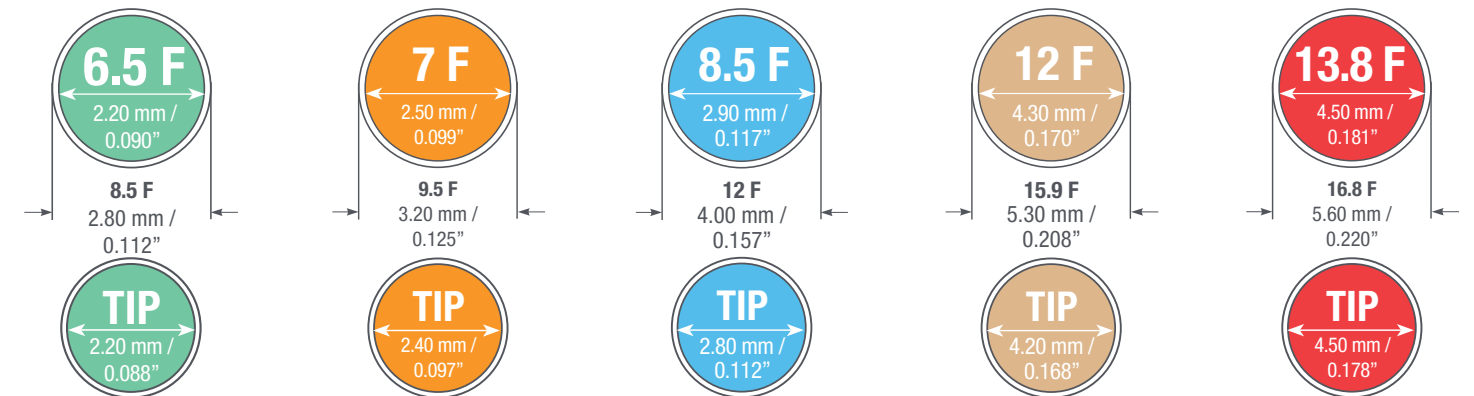

Model # Set <i>(One per box)</i>	Sheath Inner Diameter	Usable Sheath Length	Usable Dilator Length	Curve Bend
STANDARD				
DST0655507D2*	6.5 F	55 cm	73 cm	7 mm
DST0655509	6.5 F	55 cm	73 cm	9 mm
DST0659009	6.5 F	90 cm	108 cm	9 mm
DST0705509	7 F	55 cm	73 cm	9 mm
DST0705517	7 F	55 cm	73 cm	17 mm
DST0707517	7 F	75 cm	93 cm	17 mm
DST0709009	7 F	90 cm	108 cm	9 mm
DST0855517	8.5 F	55 cm	73 cm	17 mm
DST1209030	12 F	90 cm	108 cm	30 mm
DST1387139	13.8 F	71 cm	89 cm	39 mm

* Contains a second soft dilator with 80 cm usable length.

ORDERING INFORMATION

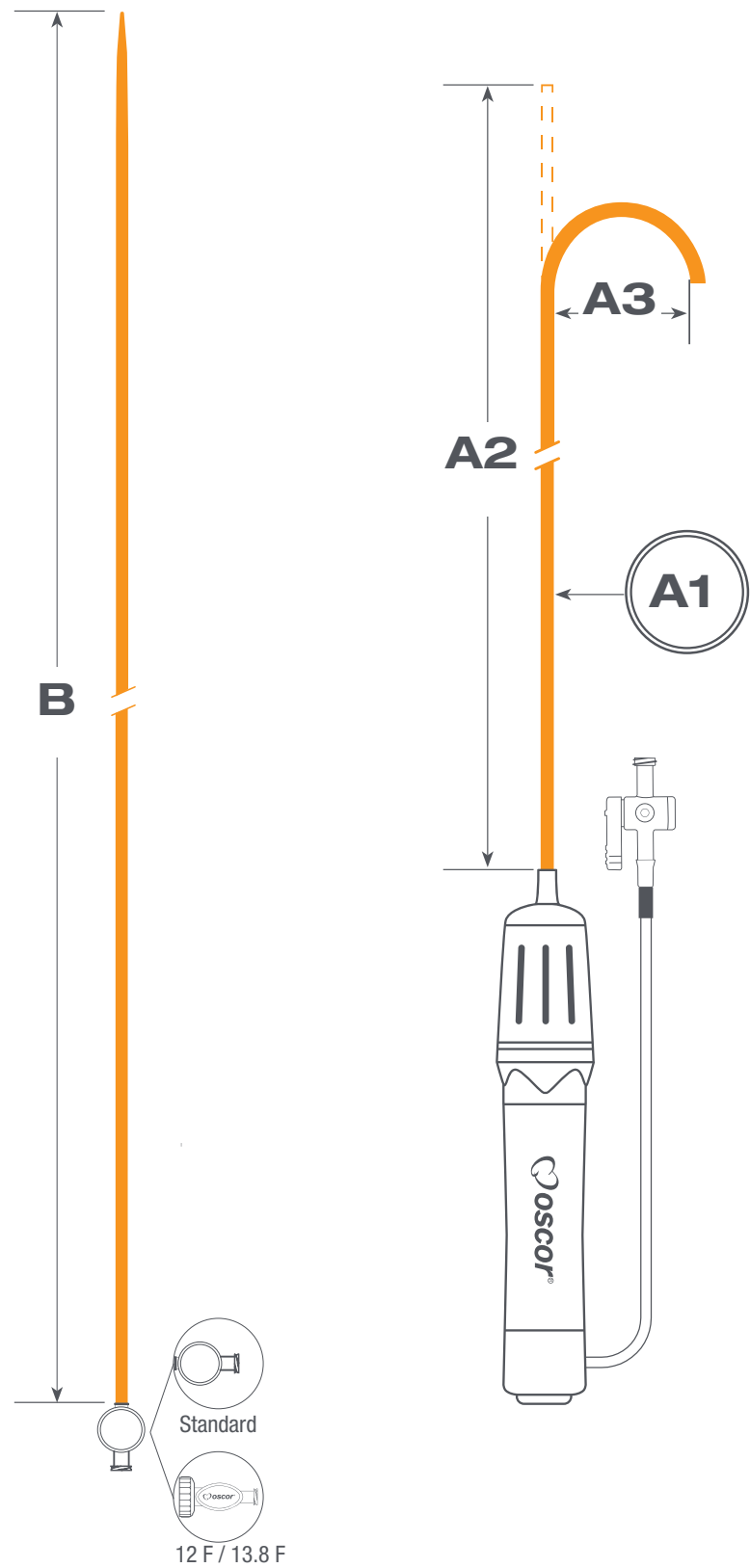

CAROTID ACCESS

Model # Set <i>(One per box)</i>	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD				
DST0659009	● 6.5 F	90 cm	108 cm	9 mm
DST0709009	● 7 F	90 cm	108 cm	9 mm



ILIAC / SFA ACCESS

Model # Set <i>(One per box)</i>	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD				
DST0654509	● 6.5 F	45 cm	63 cm	9 mm
DST0654517	● 6.5 F	45 cm	63 cm	17 mm
DST0655507D2*	● 6.5 F	55 cm	73 cm	7 mm
DST0704509	● 7 F	45 cm	63 cm	9 mm
DST0704517	● 7 F	45 cm	63 cm	17 mm

* Contains a second soft dilator with 80 cm usable length.



INFRAINGUINAL / INFRAPOPLITEAL ACCESS

Model # Set <i>(One per box)</i>	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD				
DST0654509	● 6.5 F	45 cm	63 cm	9 mm
DST0654517	● 6.5 F	45 cm	63 cm	17 mm
DST0704509	● 7 F	45 cm	63 cm	9 mm
DST0704517	● 7 F	45 cm	63 cm	17 mm



RENAL / CELIAC / SMA ACCESS

Model # Set <i>(One per box)</i>	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD				
DST0654509	● 6.5 F	45 cm	63 cm	9 mm
DST0654517	● 6.5 F	45 cm	63 cm	17 mm
DST0655509	● 6.5 F	55 cm	73 cm	9 mm
DST0655517	● 6.5 F	55 cm	73 cm	17 mm
DST0704509	● 7 F	45 cm	63 cm	9 mm
DST0704517	● 7 F	45 cm	63 cm	17 mm
DST0705509	● 7 F	55 cm	73 cm	9 mm
DST0705517	● 7 F	55 cm	73 cm	17 mm
TD				
DST0655509TD	● 6.5 F	55 cm	73 cm	9 mm
DST0705509TD	● 7 F	55 cm	73 cm	9 mm



MULTI PURPOSE ACCESS

Model # Set <i>(One per box)</i>	A1 Sheath Inner Diameter	A2 Usable Sheath Length	B Usable Dilator Length	A3 Curve Bend
STANDARD				
DST0655507D2*	● 6.5 F	55 cm	73 cm	7 mm
DST0655509	● 6.5 F	55 cm	73 cm	9 mm
DST0659009	● 6.5 F	90 cm	108 cm	9 mm
DST0705509	● 7 F	55 cm	73 cm	9 mm
DST0705517	● 7 F	55 cm	73 cm	17 mm
DST0707517	● 7 F	75 cm	93 cm	17 mm
DST0709009	● 7 F	90 cm	108 cm	9 mm
DST0855517	● 8.5 F	55 cm	73 cm	17 mm
DST1209030	● 12 F	90 cm	108 cm	30 mm
DST1387139	● 13.8 F	71 cm	89 cm	39 mm

* Contains a second soft dilator with 80 cm usable length.

Guiding sheaths are supplied sterile. For single use only.
Standard lengths are noted above. Other lengths available upon request.

Products not available in all regions.

**DISTRIBUTED BY –
BAYLIS MEDICAL TECHNOLOGIES**

2645 Matheson Blvd. East
Mississauga, ON, CA
L4W 5S4

Work: 905-948-5800 x334
Toll Free: 1-888-505-4885

baylismedtech.com



**MANUFACTURED BY –
OSCOR INC.**

3816 DeSoto Blvd.
Palm Harbor, FL 34683 U.S.A.

Phone: 727.937.2511
Fax: 727.934.9835

www.oscor.com

PRM-00875 EN J-1, 2 V-1.1. The Baylis Medical Technologies Logo is a trademark and/or registered trademark of Baylis Medical Technologies Inc. Destino is a trademark of Oscor Inc. Used under license. CAUTION: Federal Law (USA) restricts the use of the devices to or by the order of a physician. Before use, consult product labels and Instructions for Use for Indications for Use, Contraindications, Warnings, Precautions, Adverse Events and Directions for Use.

Products shown may not be approved or available in all jurisdictions.